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## Note

 <b>TECNO-GAZ</b>	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 1 of 23

# **Instruction manual for installation and use of scialytic lamp**

## **SLIM**

### **Dentist's unit**

	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 2 of 23

## Introduction

Dear installer,

we ask you to read this manual carefully before installing the product. By doing so you will be able to protect yourselves and all the possible future users against any damage.

Mark

This apparatus is a Class I medical device in accordance with the European Directive on medical devices (MDD) 93/42/EEC, May 1993, Appendix IX, and subsequent variations (including 2007/47/EC Directive, September 2007)

## Conformity

The manufacturer represents that this product complies with the essential requisites of the MDD Appendix 1 and proves such conformity by affixing the EC mark

## Manual's validity

This instruction manual for the installation of ceiling-anchor and mobile lamps is valid for the following models:

- SLIM

## Customer service

The customer service is at your disposal in case of any questions as to the device and its use, if you wish to order spare parts and for service and warranty reasons.

- TECNO-GAZ
- Strada Cavalli, 4
- I-43038 Sala Baganza - Parma - Italy
- Tél.: +39 – 0521 – 83.39.26
- Fax: +39 – 0521 – 83.33.91
- e\_mail: [info@tecnogaz.com](mailto:info@tecnogaz.com)

## Copyright

The contents of this Manual can be modified by TECNO-GAZ, without notice or further obligations, in order to include variations and improvements. Copy or translation of any part of this manual without written consent by TECNO-GAZ is forbidden.

## Modification rights reserved

TECNO-GAZ reserves the right to modify, cancel or otherwise change the data contained in this documentation at any time and for any reason

 <b>TECNO-GAZ</b>	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 3 of 23

without notice since TECNO-GAZ is constantly searching for new solutions for the development of the products. We therefore reserve the right and duty to make changes in the devices supplied in terms of form, equipment and technology.

## IMPORTANT

Please find attached a warranty coupon to be filled in and sent back to Tecno-Gaz. Such coupon completes:

- the traceability of the medical device as required by the MDD
- the Cover on the Product

the confirmation of the complete reading of the manual by the installer as a necessary and sufficient clause so that he/she is declared FIT for the installation of the equipment (having all the necessary information for the installation through this manual)


## **Manufacturer's declaration of conformity**

The Company:

**TECNO-GAZ S.P.A. Strada Cavalli n. 4 - CAP 43038 - Sala Baganza – Parma – ITALY**

declares on its own responsibility that the Medical Illumination Device for surgical and diagnostic use:

**Model:** SLIM dentist's unit      371-S2

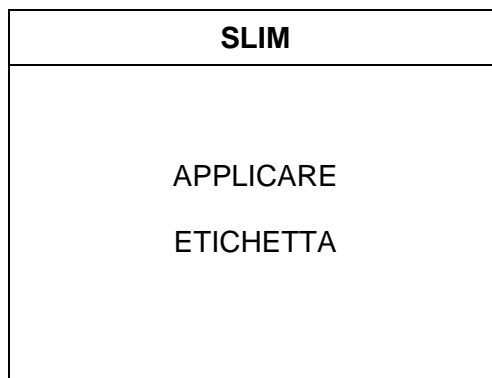constructed by TECNO-GAZ S.P.A., conforms to the Attachment VII of the 93/42/EEC Directive of May 14<sup>th</sup> 1993, and subsequent variations (including the 2007/47/EC Directive of September 5<sup>th</sup> 2007); moreover, it is in conformity with the following safety regulation:

- IEC 60601-1      (Part 1: General requirements for safety)
- IEC 60601-2-41      (Part 2: Particular requirements for safety)
- IEC 60601-1-2      (Part 1: General requirements for safety - EMC)

**Classification with reference to article 9 and Attachment IX of the 93/42/EEC Directive and 2007/47/EC**

PERIOD:      Short term (art.1, paragraph 1.1, attachment IX)

DESCRIPTION:      Non invasive Medical Device (Par.1 "Definition", art.1, paragraph 1.2, attachment IX)

Active Medical Device (art.1, paragraph 1.4, attachment IX)

CLASS:      I (Par.3 "Classification", art.1, paragraph 1.1 Rule 1, attachment IX)

Name: Paolo Bertozzi  
Position: Managing Director

Sala Baganza, 02-01-2010

## Table of Contents

<b>1</b>	<b>General information .....</b>	<b>7</b>
1.1	User's qualification .....	7
1.2	Packaging, transport, storage and characteristics of the installation site .....	7
1.3	Graphic symbols used in this installation manual .....	8
1.4	Graphic symbols used on the package .....	9
1.5	Other graphic symbols used on the device .....	9
1.6	Warranty .....	10
<b>2</b>	<b>Instructions for the device keeper .....</b>	<b>12</b>
2.1	Technical safety characteristics .....	12
2.2	Personnel instruction obligation .....	12
2.3	Warranty and responsibility .....	12
2.4	Structural changes or variations.....	12
2.5	Disposal after use.....	13
<b>3</b>	<b>Product installation .....</b>	<b>13</b>
3.1	Lamp installation .....	13
3.2	First switch on .....	13
3.3	Installation check and scialytic lamp testing operations before use .....	14
<b>4</b>	<b>Importance of personal safety.....</b>	<b>15</b>
4.1	Intended use.....	15
4.2	Environmental conditions.....	15
4.3	Technical safety conditions.....	15
4.4	Other safety conditions (secondary effects).....	16

<b>5</b>	<b>Lamp description and operation .....</b>	<b>16</b>
5.1	Description of the device .....	16
5.2	Description of the operation .....	17
<b>6</b>	<b>Cleaning and disinfection .....</b>	<b>17</b>
6.1	Cleaning the device.....	17
6.2	Disinfection.....	18
6.3	Sterilizing the handles .....	19
6.4	Yearly inspections by the keeper .....	20
<b>7</b>	<b>Adjustments .....</b>	<b>20</b>
7.1	Adjustment of the swing arm.....	20
7.2	Adjustment of the braking force .....	20
7.3	Search for failure .....	21
7.4	Spare part list .....	21
<b>8</b>	<b>Technical data.....</b>	<b>21</b>
	<b>Notes.....</b>	<b>23</b>

 <b>TECNO-GAZ</b>	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 7 of 23

## 1 General information

### 1.1 User's qualification

Personnel's qualification The following scheme describes in chronological order the life of the scialytic lamp and the qualifications required by TECNO-GAZ:

Installation	TECNO-GAZ or Qualified Installer
Use	Qualified/authorized medical personnel
Ordinary maintenance	TECNO-GAZ or maintenance by the qualified technical body
Extraordinary maint.	TECNO-GAZ or Qualified Maintenance Technician
Service	TECNO-GAZ or Authorized Retailer
Cleaning	Duly trained personnel
Demolition	TECNO-GAZ, Qualified Installer or Maintenance Technician

#### **DANGER**

**TECNO-GAZ refuses all responsibility for any damage to people or things deriving from the use of the product by non-qualified operators.**

**By qualified Installer/Maintenance Technician we mean the person/s who has/have attended a course at TECNO-GAZ or read this manual carefully and, as confirmation thereof, the warranty coupon attached must have been sent back**

### 1.2 Packaging, transport, storage and characteristics of the installation site

Packaging	1 box containing a complete scialytic lamp with installation instruction manual and user's manual.
Transport	Transport is carried out by any road haulage contractor as long as they respect the following characteristics: Temperature (°C): -15 / +60 Humidity: 10 / 75 % Atmospheric pressure (h/Pa): 500 / 1060
Storage	The devices packaged must be stored (warehoused) in a dry place and at the following temperature: Temperature (°C): -15 / +60 Humidity: 10 / 75 % Atmospheric pressure (h/Pa): 500 / 1060

	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 8 of 23

#### Installation site

The site appointed for the installation of the equipment must have the following characteristics:

Temperature (°C): +10 / +40

Humidity: 30 / 75 %

Atmospheric pressure (h/Pa): 700 / 1060

### 1.3 Graphic symbols used in this installation manual

In these installation instructions and on the lamp important directions are marked by means of symbols and identifying words.

Identifying words such as **DANGER**, **WARNING** or **ATTENTION** indicate the classification of the risk of damage. This is highlighted with various symbols.

**DANGER** signals an immediately dangerous situation which might cause death or serious damage.

**WARNING** signals a potentially dangerous situation which might cause death or serious damage.

**ATTENTION** signals a potentially dangerous situation which might cause moderate or slight damage.

The following triangular symbol coupled with a side explanation indicates which danger is being faced:

**Electric shock, Mechanical danger due to hanging masses (quick release of a suspension arm)**

## 1.4 Graphic symbols used on the package

The symbols present on the boxes of the packages are listed below:

High side of the package

Maximum number of packages stackable

Breakable package

Humidity-suffering package

Do not overlap packages with pallet

Weight of the package

## 1.5 Other graphic symbols used on the device

The symbols present on the lamp are listed below:

B-Type device. Indicates the level of protection against direct and indirect contact

Graphic symbol proving the EC marking of the product

Symbol indicating the manufacture date (month and year)

Fuses used by the device

	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 10 of 23

## 1.6 Warranty

Cover on the product	1. The product is covered by guarantee for a period of 12 months, including electric parts.
Starting date of the warranty	2. The warranty becomes effective from the lamp installation date at the customer's premises exclusively if the installation is carried out by Tecno-Gaz personnel. In any other case the warranty becomes effective from the lamp dispatch date from Tecno-Gaz warehouse to the customer's premises.
Transport Document Date	3. In case of contention, the date specified on the "Transport Document" given with the goods is considered valid.
Repairs	4. Repair or replacement, under guarantee, of a part is carried out for well-assessed manufacture reasons and at the unchallengeable discretion of Tecno-Gaz. The warranty does not include personnel's travels, packaging and transport charges.
Standard wear	5. The warranty does not include components subject to normal wear and tear (including but not limited to: halogen lamps, fuses, relays, bearings, etc.)
Exclusion	6. The warranty does not include: <ul style="list-style-type: none"> <li>- failure and faults due to improper installation;</li> <li>- failure or defects caused by carelessness, negligence, abnormal use of the device or other causes not attributable to the constructor;</li> <li>- interventions due to presumed defects or convenient inspections.</li> </ul>
Replacement	7. The right to replace the whole Lamp is not acknowledged. 8. The warranty does not grant any indemnity for direct or indirect damage of any nature to people or things, due to the inefficiency of the lamp.
Improper use	9. Tecno-Gaz can not be held responsible for failure or damage caused by improper use of the product, non-maintenance or neglect of the basic principles of good maintenance (negligence).
Use stop	10. No indemnity is acknowledged for Lamp stop.
Loss of warranty	11. The warranty is lost automatically if the lamp is tampered with, repaired or modified by the buyer or third parties not authorized by Tecno-Gaz.

 <b>TECNO-GAZ</b>	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 11 of 23

Technical interventions	12. For the interventions, the buyer must apply to the retailer only or the service centers as recommended by Tecno-Gaz.
Replacing components	13. The components replaced under guarantee must be returned to Tecno-Gaz carriage free.
Non-return	14. The non-return implies the charge of the cost of the part to the applicant.
Returns	15. Tecno-Gaz does not accept returns by final users.
Repairs	16. The returns to Tecno-Gaz for repair must be managed by the retailer or the service center chosen by the final user according to the commercial procedure CM-P-003 called "MANAGEMENT OF REQUEST FOR RESTITUTION, REPAIRING, REPLACEMENT OF TECNO-GAZ EQUIPMENT AND CLIENT COMPLAINTS".
Authorization for returns	17. The returns to TECNO-GAZ must be documented and authorized as per internal procedures.
Technical documents	18. The products returning to Tecno-Gaz must be attached with the return authorization documents and a document describing the failure.
Return dispatch	19. All products under repair must be sent to Tecno-Gaz carriage free and suitably packaged (it is compulsory to use the original package).
Non-responsibility	<p>20. The constructor does not deem himself responsible for the safety, reliability and performance of the device if:</p> <ul style="list-style-type: none"> <li>- the installation, additions, settings, repairs are not carried out by people authorized in advance by Tecno-Gaz;</li> <li>- the electrical system of the installation site does not comply with the national laws and/or regulations in force for electrical systems in sites used as medical premises and similar norms;</li> </ul> <p>the equipment is not used in accordance with the instructions for use.</p>

 <b>TECNO-GAZ</b>	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 12 of 23

## 2 Instructions for the device keeper

### 2.1 Technical safety characteristics

Cleaning personnel	The cleaning and disinfection operations of the scialytic lamp system described in Chapter 6 must be carried out exclusively by duly instructed cleaning personnel.
Service personnel	The control and maintenance operations described in Chapter 7 must be carried out exclusively by qualified service personnel.

### 2.2 Personnel instruction obligation

Instructing the users	<p>Instruct the personnel following the instructions for use as regards the lamp control, cleaning and care operations.</p> <p>Such personnel must be instructed by the keeper according to what is specified in this manual.</p>
-----------------------	-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

### 2.3 Warranty and responsibility

TECNO-GAZ refuses all responsibility for the unreliable working of the scialytic lamps in case:

- the installation, the changes and the repairs are not carried out by a TECNO-GAZ service technician or personnel authorized by TECNO-GAZ,
- the device is not used in accordance to its intended use, the norms and the instructions for use.

### 2.4 Structural changes or variations

Arbitrary changes	For safety reasons, all structural changes or variations made in the device on one's own initiative are forbidden. In case of similar changes or transformations, the warranty of the producer for this device is annulled. This way the producer does not assure any warranty for damage or injuries caused by structural changes or variations made in the device on one's own initiative or by the use of non-original spare parts.
Use TECNO-GAZ original spare parts only	Using spare parts not provided by TECNO-GAZ or its representatives implies annulment of the warranty.

	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 13 of 23

## 2.5 Disposal after use

Life end disposal

The used device contains valuable materials which can be reused.

Dispose of the used device respecting the environment and in accordance with the national waste disposal directives in force.

## 3 Product installation

**Before starting the installation make sure that no package is missing, in bad condition or damaged during transport and that the contents conform to what is specified above.**

**Complaints will be considered only if the seller or the forwarder are warned immediately. Each complaint must be made in writing.**

**Goods are always carried at the buyer's risk.**

**Keep the original package in case the lamp is to be returned.**

### 3.1 Lamp installation

**ATTENTION**

Incorrect levelling of the stem may cause the horizontal arm of the lamp to make an undesired spontaneous move, when passing from the non-balance condition to the balance one.

See drawing 76

Align the horizontal arm pin (1) to the stem of the dental (2).


Insert the electrical cables (3), coming from the lamp structure, along the dental stem (2), making them coming out from the base.

Insert completely the pin in its seat.

Connect the power supply (4) to the red and black cables coming from the structure, through connectors (5) and (6), be sure to maintain the connect polarity.

### 3.2 First switch on

Before supplying the lamp with voltage, remove the protection film from the head screen surface and from the floor base cover.

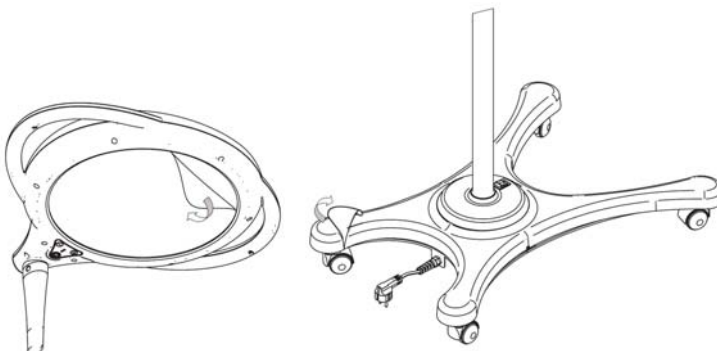	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 14 of 23

## ATTENTION

Now it is possible to power the lamp to make sure it gives out light. To do so you must:

- 1 – close the system circuit by means of the head system switch (customer's responsibility)
- 2 – press the green switch on the base cover (floor version)
- 3 – press the I/O keyboard

## 3.3 Installation check and scialytic lamp testing operations before use

## ATTENTION

The following notes are to be considered compulsory during the installation check as they prove the correct verification of all the points set out. For such reason every single point must be ticked when it is treated.

1. Make sure the dental stem is suitable for installing the lamp ف
2. By means of a spirit-level, make sure the dental stem is perpendicular to the ground ف
3. Check the device earthing and make sure the earthing terminals are well tightened ف
4. Make sure the movement mechanism works perfectly. Check the mechanical operation of the lamp by directing and rotating it and make sure that all movements are frictioned so that the lamp is stable in any position and does not lose its position ف
5. After turning on the lamp, it must give out light from the head ف

Seal and signature of the installer:

---

 <b>TECNO-GAZ</b>	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 15 of 23

## 4 Importance of personal safety

### 4.1 Intended use

Use complying with the norms

The SLIM scialytic lamp is used to light the area of the patient undergoing a surgical operation or an examination and intended to be used in operating rooms or medical clinics.

Work field

The work field is a distance from about 70 to 140 cm from the operation point.

#### Single lamp (minor surgical luminaire):

Definition

In accordance with EN60601-2-41, a SLIM 12 single lamp is a minor surgical luminaire and can be used only on operations where the interruption of the lighting does not cause any risks for the patient.

### 4.2 Environmental conditions

- The SLIM 12 scialytic lamp system is not suitable for the use in explosion-risky areas.
- The SLIM scialytic lamp system is not suitable for the use in the presence of inflammable mixtures of anesthetics with air, oxygen or NO<sub>2</sub> (laughing gas).
- During working, the room temperature must be between 10°C and 40°C.
- The relevant humidity must be between 30% and 75%.
- The atmospheric pressure must be between 700 and 1060hPa.

### 4.3 Technical safety conditions

The safe use and proper working of the SLIM lamp is guaranteed if:

Safe anchorage anchorage of the lamp to the dental is safe from the static point of view and there is a static stability test,

Systems the electrical systems of the sites concerned comply with the local norms in force,

Authorized personnel changes or maintenance of the lamp are made by TECNO-GAZ authorized personnel,

	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 16 of 23

Proper installation and commissioning	the lamp has been installed following the currently valid installation instructions and has been commissioned regularly by the authorized installer,
Original spare parts	TECNO-GAZ original parts only are used for service, repair and structural variation operations and as accessories.

#### 4.4 Other safety conditions (secondary effects)

Optical security	Do not direct the luminous source in the eyes of patient or operator, in order to avoid possible dazzling effect and eventual retinal damages. PRECAUTIONS: - Do not direct the luminous source in patient's or operator's eyes - It is compulsory to cover the patient's eyes with an adequate protection
Improper use	Do not lay any object on the body of the lamp or hang objects on the arm or the body of the lamp, since a safe positioning is no longer ensured and such objects may fall on the operation area. The load of heavy objects (body weight of a human being) can damage the mechanics.
Cover the heads	During working the head of the lamp must not be covered since it would prevent the thermal exchange with the environment and cause overheating of the lamp.
Collisions	Collisions of the mobile arms and the lamp body must be avoided. In case of violent collision, the lamps may be damaged or pieces of paint may chip and fall on the operating field in the patient area.

## 5 Lamp description and operation


### 5.1 Description of the device

Versions	This version of the SLIM lamp is related to an: - Dentist's unit
See drawing 77	DENTIST'S UNIT version: anchoring pin (1), horizontal arm (2), heeling arm (3), lamp head (4), function control keyboard (5), sterilizable handle (6), power supply (7).

## 5.2 Description of the operation

### Control panel

The function control keyboard is applied to the scialytic lamp head. Such keyboard allows to turn the lamp on and off by means of the membrane key I/O (1), to regulate the light intensity by pressing the keys with the symbols '+' (2) and '-' (3), and display the level of intensity reached through 3 green position micro-LEDs (4).
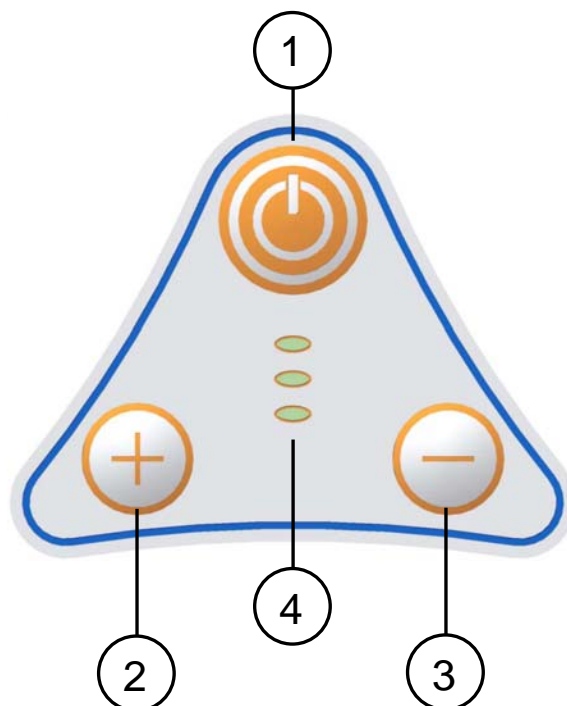
### Light field

The light field is fixed and therefore any adjustment of its diameter dimension is prevented.

## 6 Cleaning and disinfection

### 6.1 Cleaning the device

	<b>WARNING – Electric shock danger</b>
-------------------------------------------------------------------------------------	----------------------------------------

Turn off the lamp by the master switch of the operating room and make sure it does not turn on again; in the case of floor versions remove the power supply plug from the electric system socket, protect the device against water sprinkling and do not clean/disinfect it with liquids.

Clean with alcohol, low-alkaline and chlorine-less detergent.

 <b>TECNO-GAZ</b>	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 18 of 23

## ATTENTION

Do not use abrasive products, petrol, paint thinners, alkaline detergents, acids;  
proportion the detergents so that so that no liquids seep into the lamp bodies and the support arm system.  
Clean the lamp with a damp but not wet cloth.

## 6.2 Disinfection

	<b>WARNING – Electric shock danger</b>
-----------------------------------------------------------------------------------	----------------------------------------

Turn off the lamp by the master switch of the operating room and make sure it does not turn on again; in the case of floor versions remove the power supply plug from the electric system socket,  
protect the device against water sprinkling and do not clean/disinfect it with liquids.  
Let the lamp body cool down. Disinfect the lamp body only when it is cold.  
Use alcohol or aldehyde-based products.

## WARNING

Disinfectants may contain substances harmful to the health: use disinfectants only respecting the health standards fixed by the hospital, the device keeper must respect the requisites imposed by the national competent hygiene and disinfection commission.

## ATTENTION

To avoid damage to stainless steel, aluminum parts, only use disinfectants with no chlorine or halogens,  
to prevent plastic parts from becoming fragile, only use low-alcohol disinfectants,  
proportion the detergents so that so that no liquids seep into the lamp bodies and the support arm system.  
Clean the lamp with a damp but not wet cloth.

### 6.3 Sterilizing the handles

	<b>WARNING – Danger for the patient</b>
-----------------------------------------------------------------------------------	-----------------------------------------

Replace the handles as soon as they show cracks or deformations, since they might fall on the area of the wound.

The lamp keeper must respect the requisites imposed by the national competent hygiene and disinfection commission.

Installation / uninstallation of the handle:

- press the two lateral safety keys of the handle and extract the handle itself.
- insert the handle till the bottom until the safety keys don't fall in their seats.

Cleaning, disinfection and sterilization of the handle

the handles are made of plastic material resistant to heat and collisions (PSU).

They can be cleaned with a detergent having medium alkaline level but no active chlorine.

To disinfect the handles we recommend to use alcohol or aldehyde-based products. The disinfectants must be approved by the manufacturer for use on the polysulfone (PSU).

Before sterilizing, rinse the handles.

The handles can bear about 300 steam sterilization cycles if the following instructions are respected:

- steam sterilization at a121°C 1,3bar from 25 to 30 minutes,
- or

- steam sterilization at a134°C 2,3bar 4 minutes.

Put the handles upright with downward open side.

Do not exceed the sterilization temperature of 134°C.

Avoid contact of the handles with other objects when disinfecting.

**Each device, in time, is subject to some wear and tear. The safety and working of the device must then be checked in the inspection and maintenance intervals.**

## 6.4 Yearly inspections by the keeper

	<b>ATTENTION</b>
-----------------------------------------------------------------------------------	------------------

Respect the yearly inspection intervals and check the device as follows:

- damage to the painting with danger in case of open wounds
- cracks on the plastic part
- detachment of parts
- integrity of the earthing circuit

## 7 Adjustments

### 7.1 Adjustment of the swing arm

See drawing 65

The lamp is supplied already balanced and needs no further adjustment. If, in time, the spring balancing swing arm stiffens or loosens, it is possible to intervene mechanically adjusting the compression of the internal spring. Slide the rubber grommet forward (1) and the cover (2) along the heeling arm (3). Insert a pin (4) having maximum diameter of 4mm into the holes of the ring nut (5) and rotate in the directions of the arrows to increase/decrease the charge of the spring.

If the heeling arm lowers, the elastic force of the spring is insufficient:

- rotate the lever downwards to charge the spring.

If the heeling arm still goes up, the elastic force of the spring is too high:

- rotate the lever upwards to discharge the spring.

At the end of the adjustment, put the cover back to its original position.

### 7.2 Adjustment of the braking force

See drawing 65

As for all the mechanic parts, brakes also are subject to wear and tear.

If the lamp body does not automatically keep the position it is put to, it is necessary to adjust the braking force by acting on the screws of the brakes.

Horizontal arm brakes

Use a 2.5 hexagonal wrench (7) to increase the braking force, rotating the arm brake dowels (6) clockwise.

### 7.3 Search for failure

n	Problem	Solution
1	The lamp does not keep the position	Check the leveling of the dental stem Tighten the brakes to the articulations in order to increase the frictions.
2	The lamp does not work	Check the the hooking of the electrical connectors. Check the presence of voltage in the lamp head
4	The lamp flickers	Contact the service center
5	The light beam on the operating field is not on focus (faulty meeting of the light fields through handle rotation)	Contact the service center.
6	The lamp does not turn on.	Check the power supply voltage. Electronics are damaged: contact the service center

### 7.4 Spare part list

Description	Order code
Sterilizable handle	Z180045
Electronic card SLIM	Z300617
Membrane keyboard	Z300220

## 8 Technical data

Technical data on light	SLIM
Illumination $E_c$ at a distance of 1 m [Lux] $\pm 10\%$	50.000
Color temperature [K] $\pm 5\%$	5.000 K ( $\pm 5\%$ )
Color rendering index $R_a$ [-]	90

Diameter of the light field d <sub>50</sub> [mm]	85
Diameter of the light field d <sub>10</sub> [mm]	140
Illumination depth [mm]	1500
Maximum irradiation [W/m <sup>2</sup> ]	119
Irradiation / Illumination [mW/m <sup>2</sup> lx]	2,38
Maximum irradiation in the UV [W/m <sup>2</sup> ]	0,001
Focus by handle	No
<b>Data on electrical connection</b>	
Primary alternating voltage [Volt ac]	100÷240
Frequency [Hz]	50/60
Absorbed power [VA]	40
Light source	N°12 Led x 2W
LED light source duration [h] (this datum can vary according to voltage peaks and the frequency of use)	50.000
Light intensity control	50 - 100
<b>General data</b>	
Color	RAL 9003
Directive	2007/47/EC
Standard	EN 60601-2-41
Electrical safety class	Class I
Protection against direct and indirect contacts	B-type device
<b>Dimensions</b>	
Lamp body diameter [cm]	30
Diameter of the poly-elliptical reflector [cm]	5,5

 <b>TECNO-GAZ</b>	<b>Instruction manual for use and installation</b>	MIU_011_GB	24/03/10
		Rev.0	Page 23 of 23


Light emission surface [cm <sup>2</sup> ]	285
Scialytic dental lamp weight [Kg]	10
<b>Certificates</b>	
	Complying with directive 2007/47/EC and 93/42/EEC


## Notes


CEILING PLATE SHEET TEMPLATE

(SCALE 1:1)


*DIMA PIASTRA A SOFFITTO*

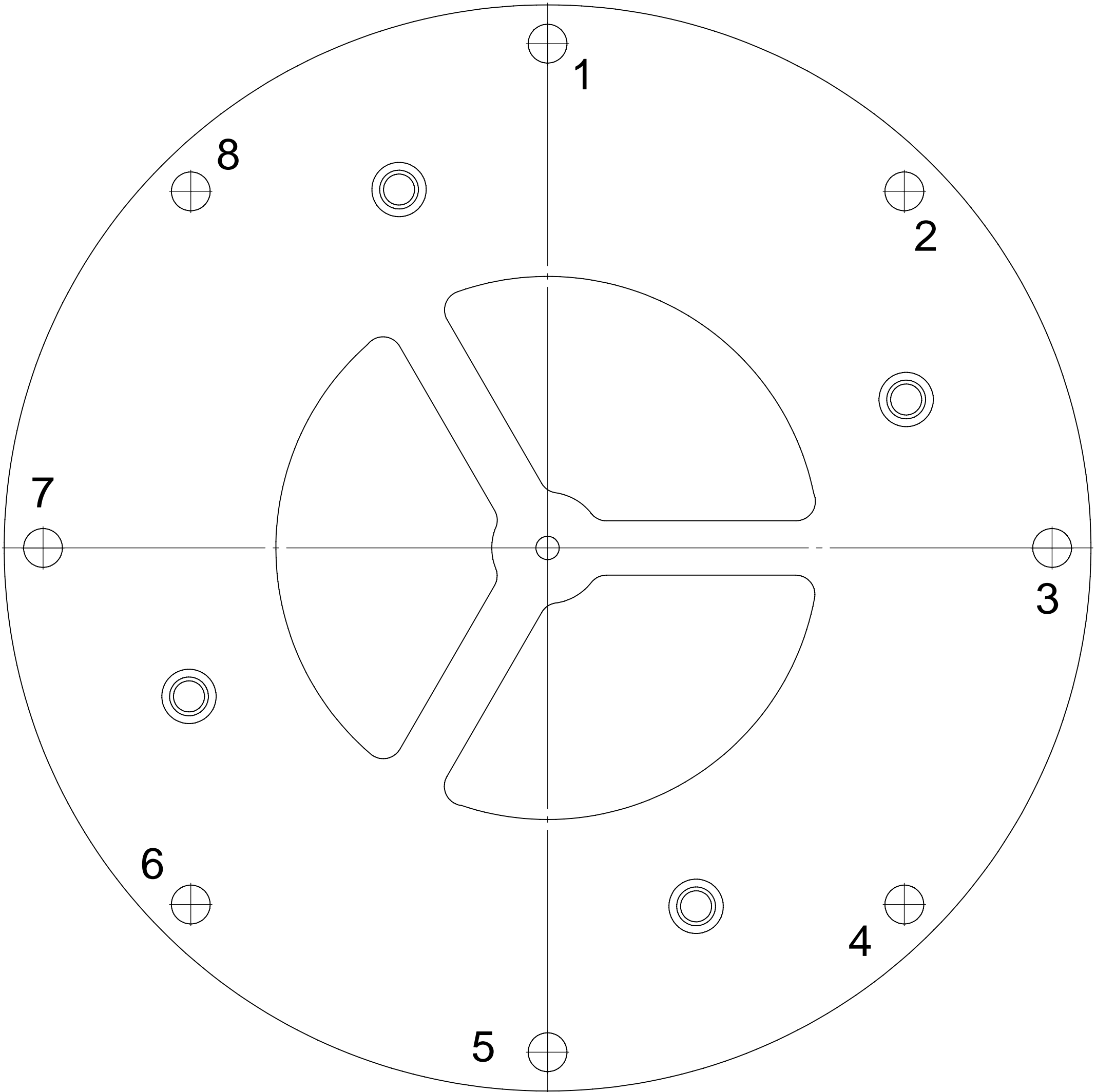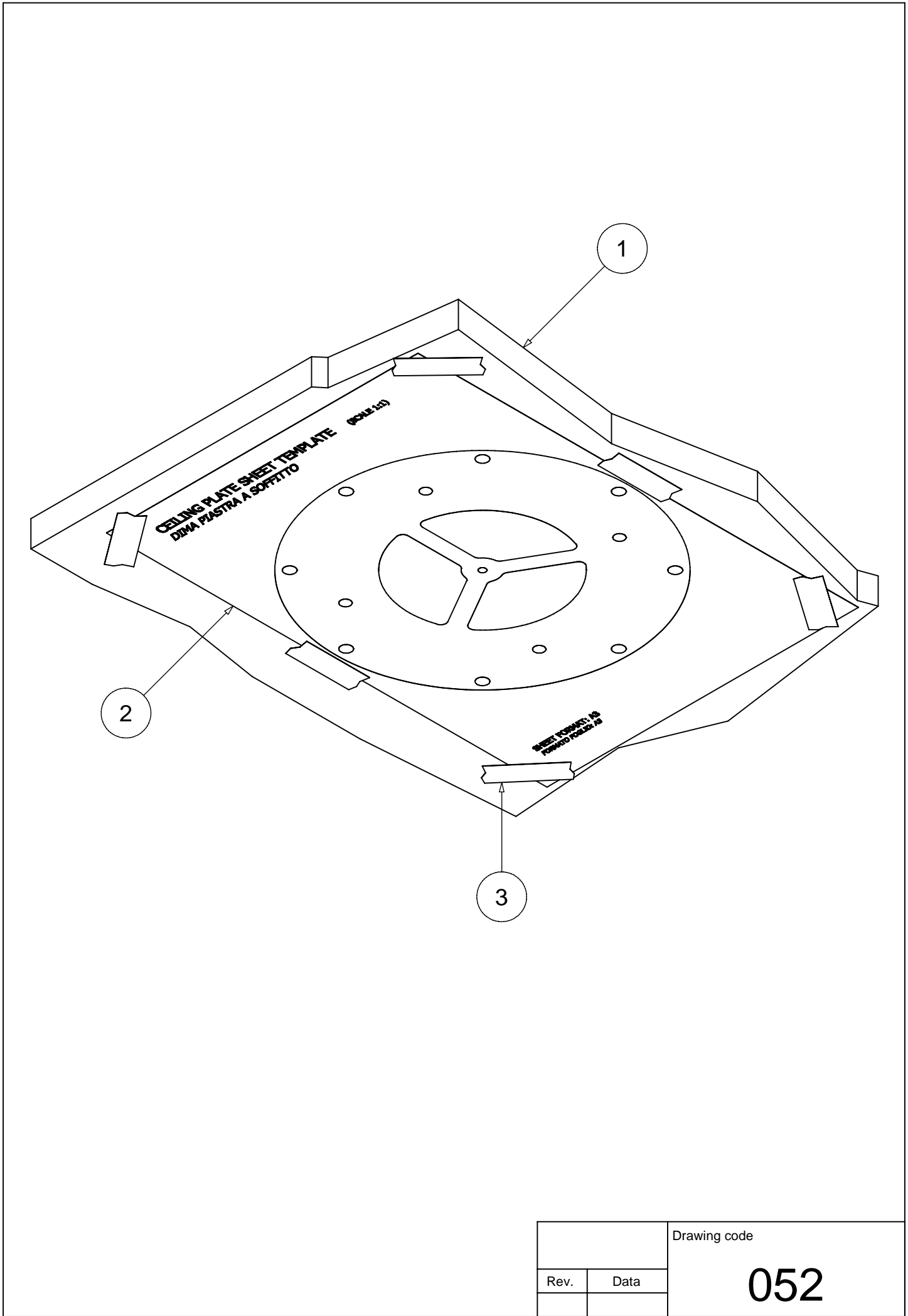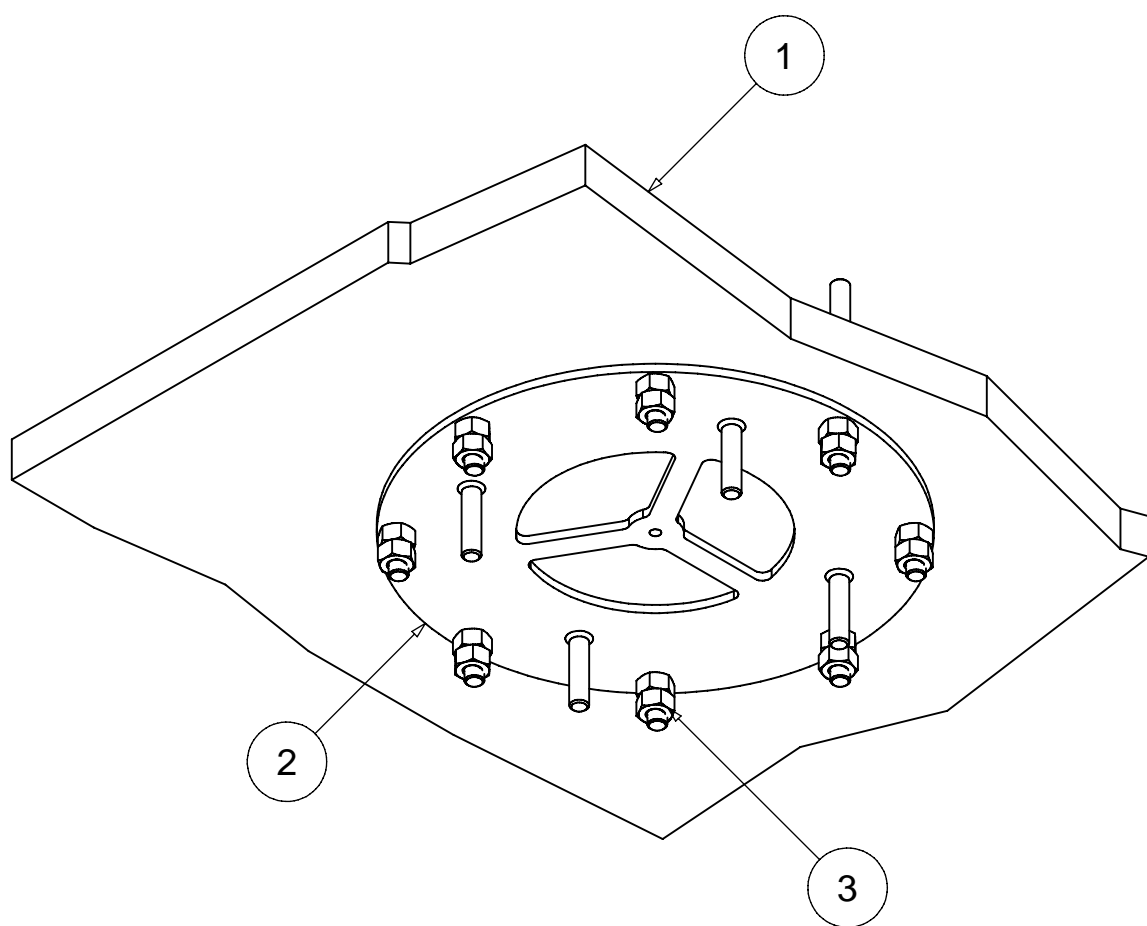
		Drawing code  <b>053</b>
Rev.	Data	

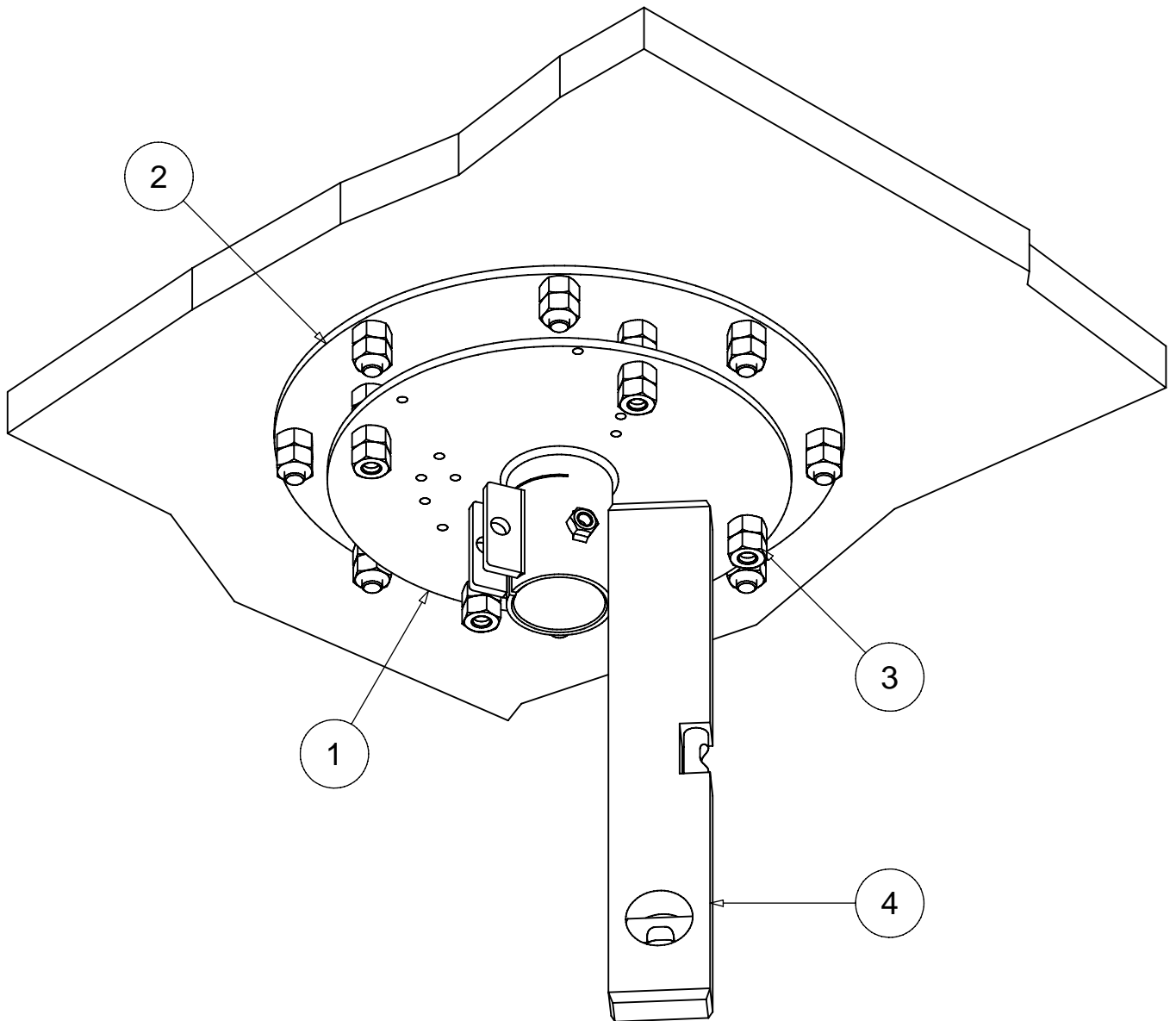
		Drawing code  <b>054</b>
Rev.	Data	

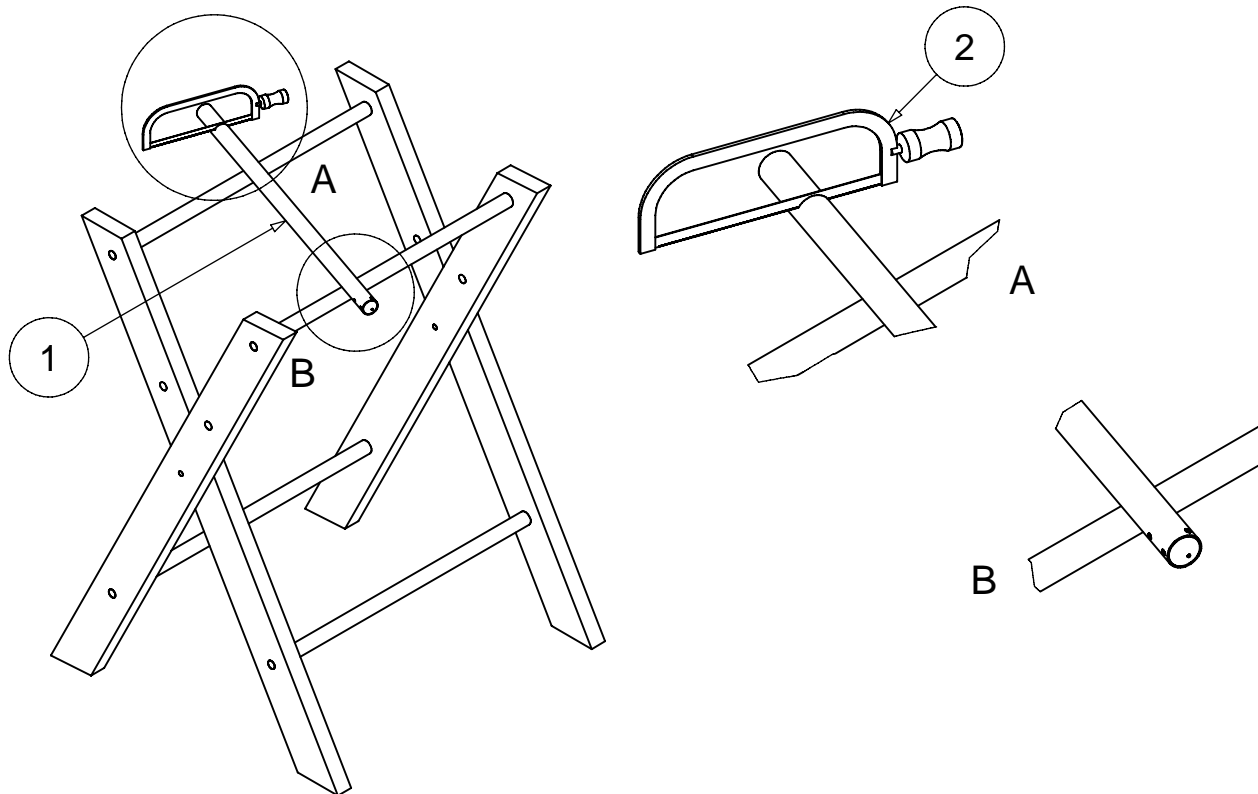
## Tube length list [mm]

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

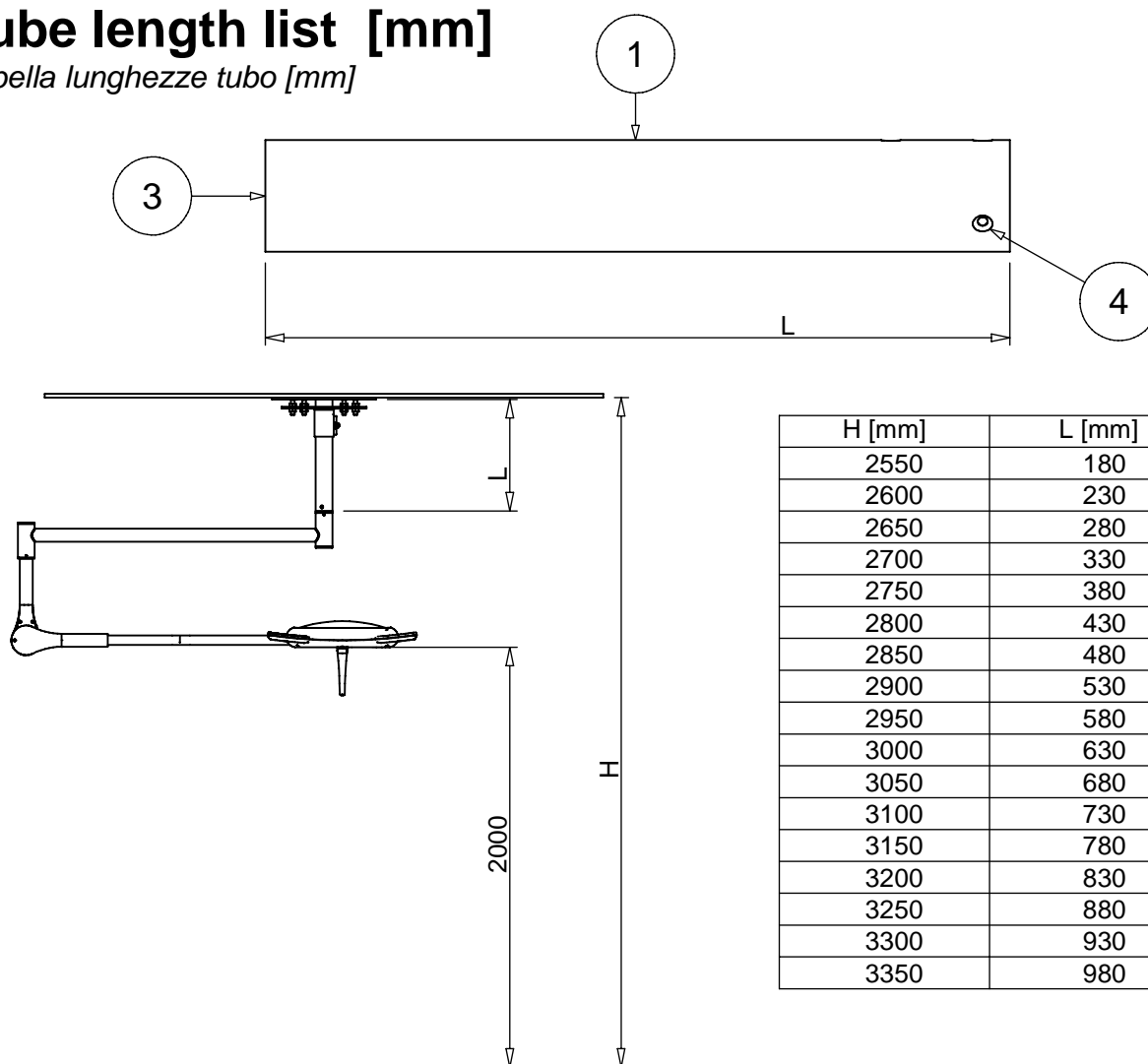
		Drawing code
Rev.	Data	055

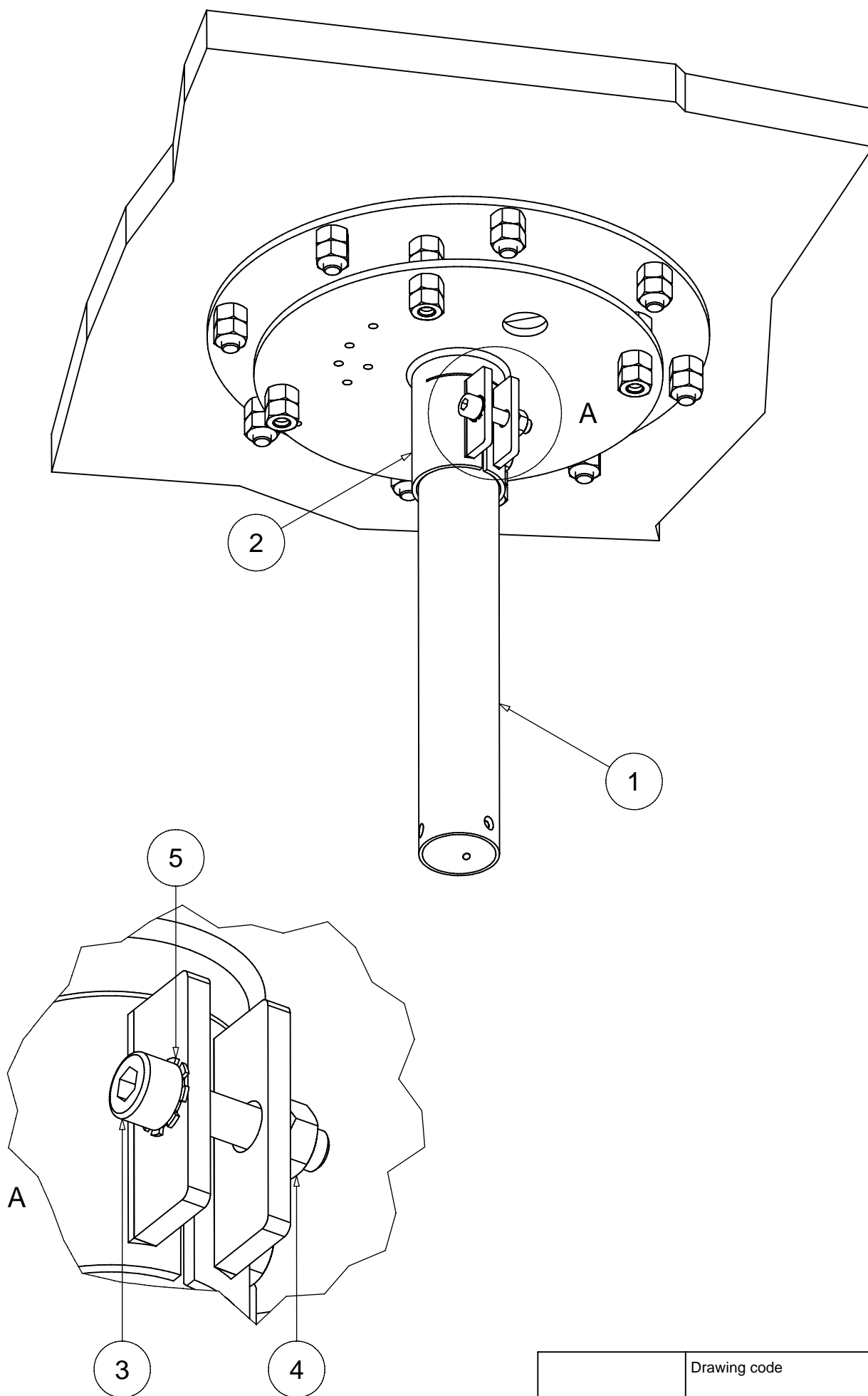
		Drawing code <b>056</b>
Rev.	Data	

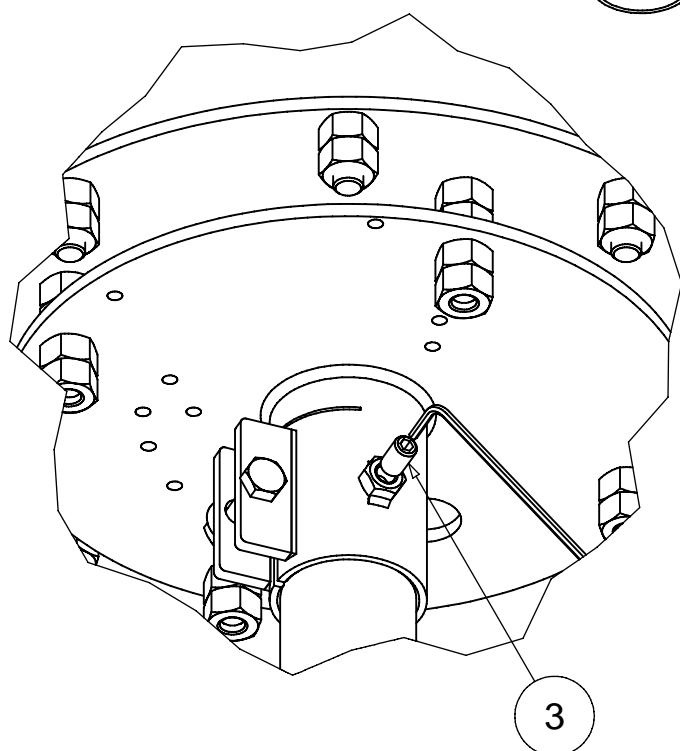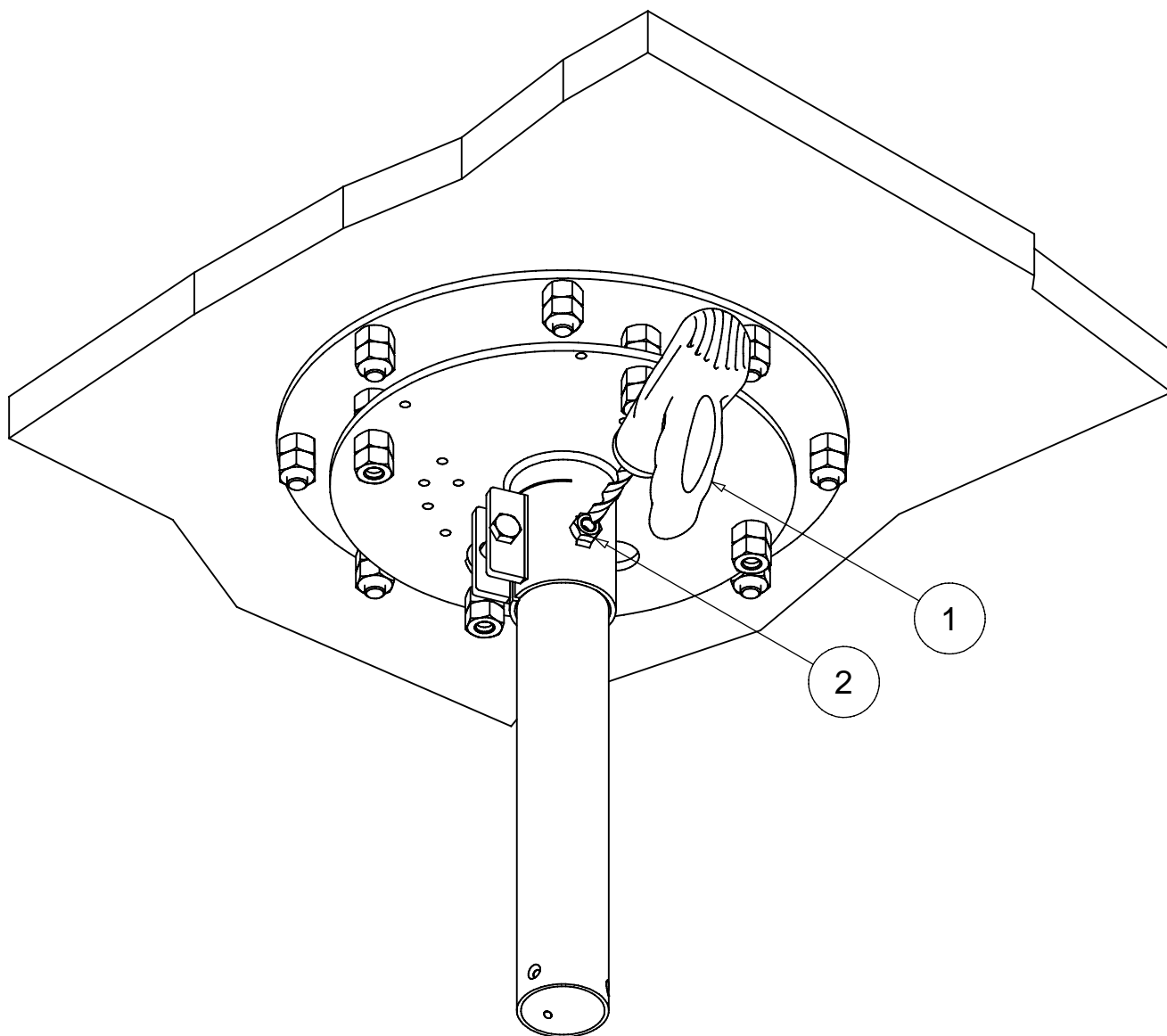
		Drawing code <b>057</b>
Rev.	Data	




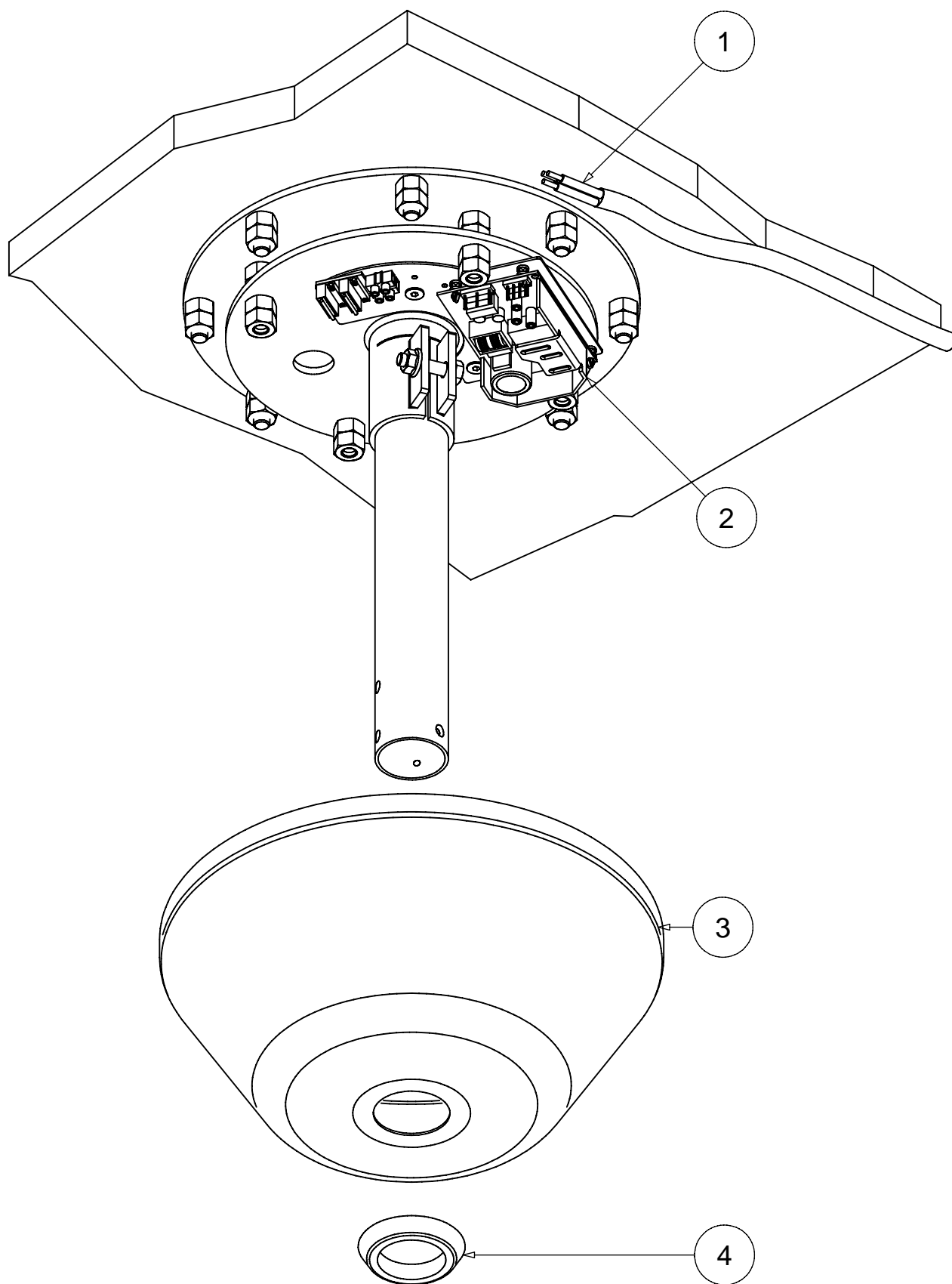
		Drawing code <b>058</b>
Rev.	Data	

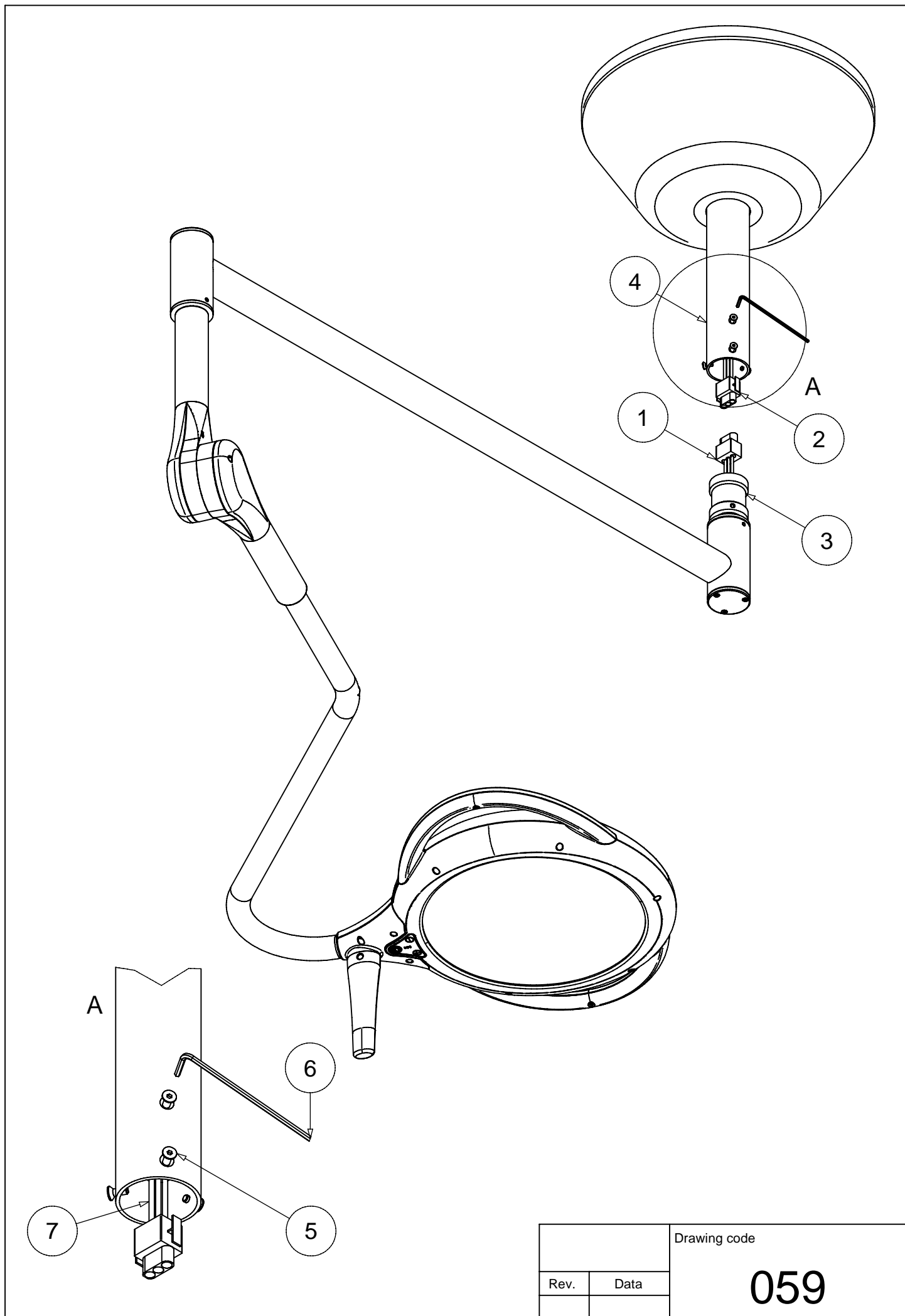
		Drawing code <b>059</b>
Rev.	Data	

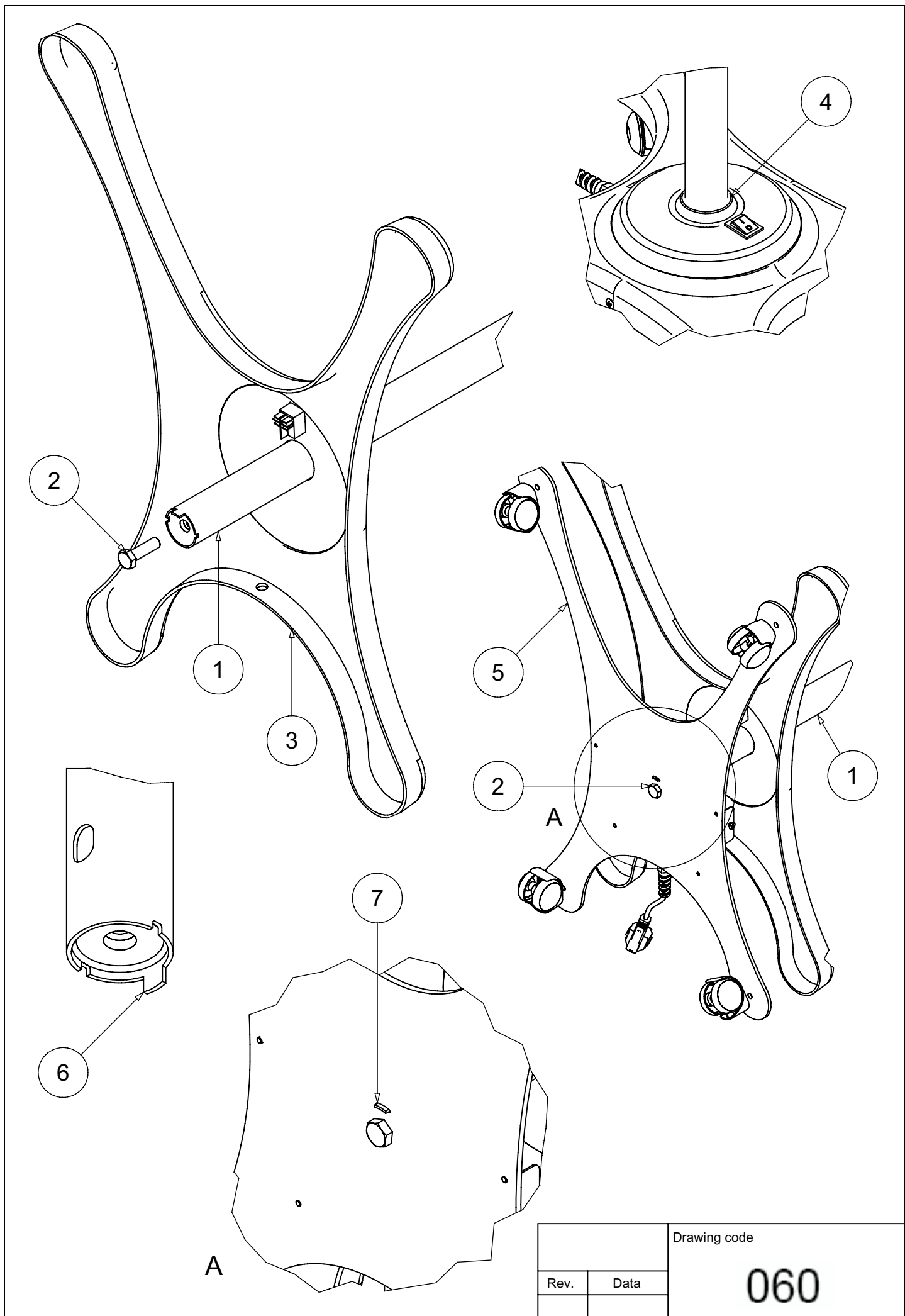
		Drawing code
Rev.	Data	060

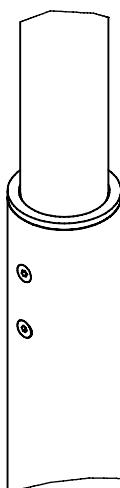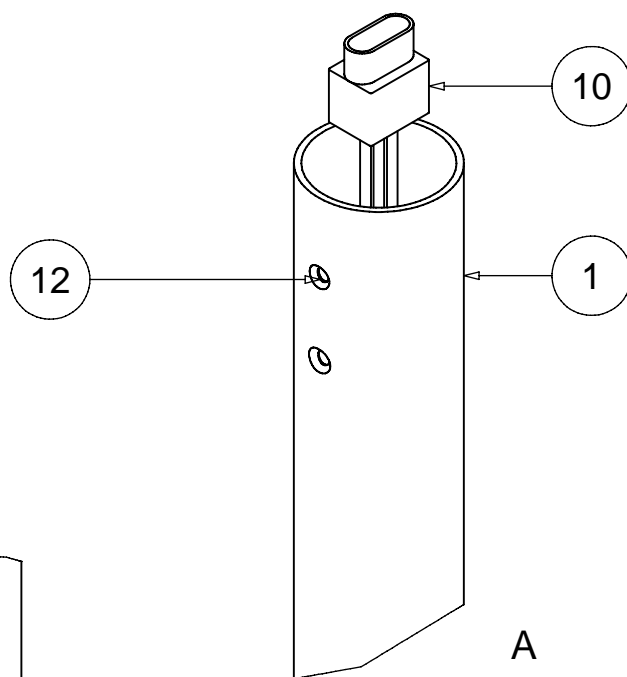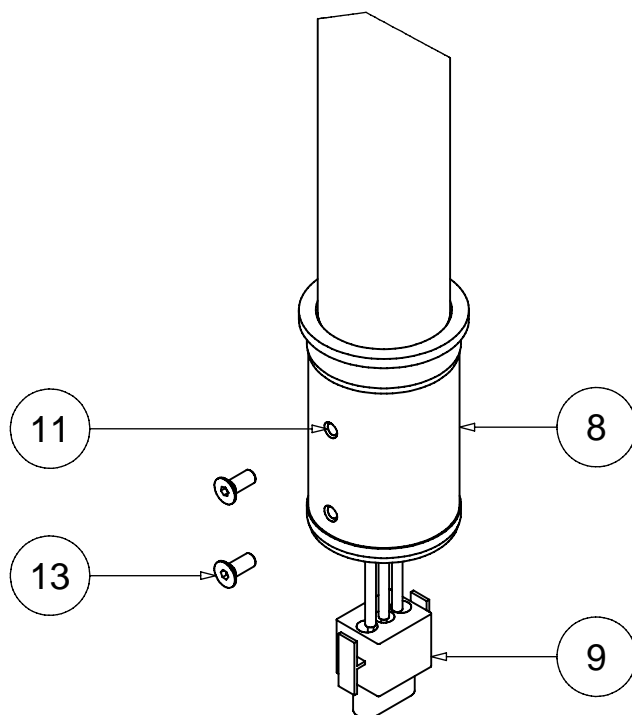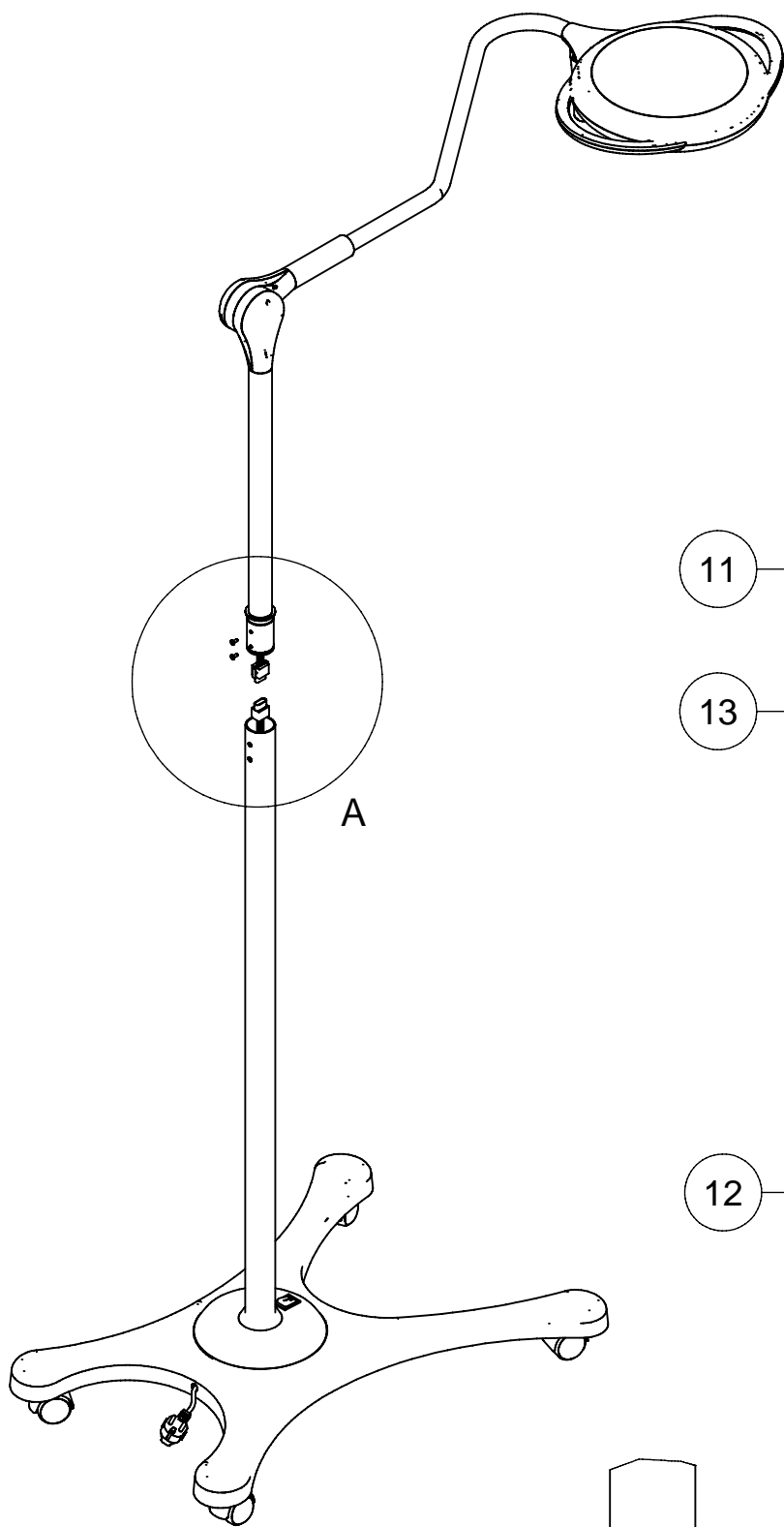
		Drawing code <b>061</b>
Rev.	Data	

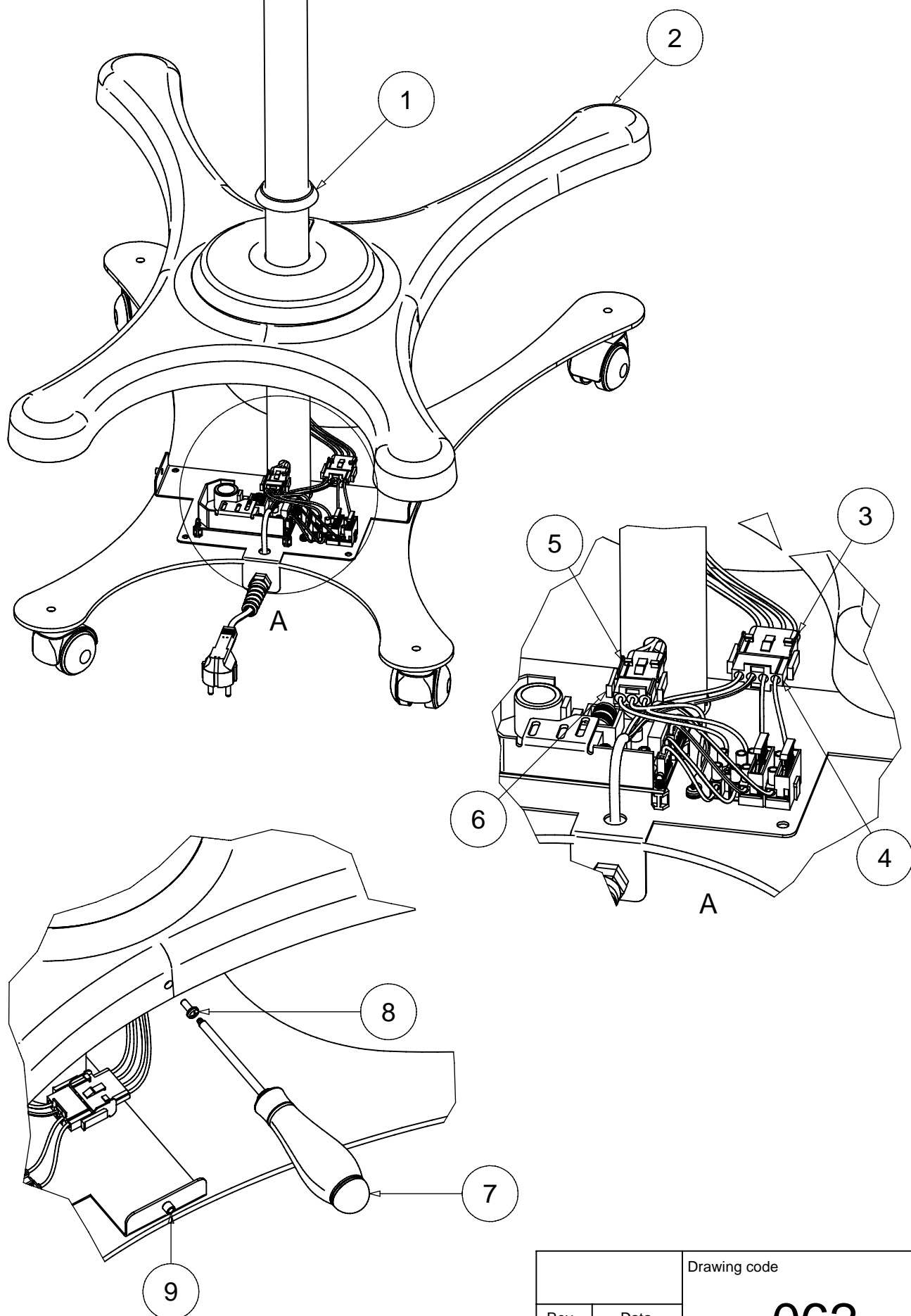
		Drawing code	
Rev.	Data	062	

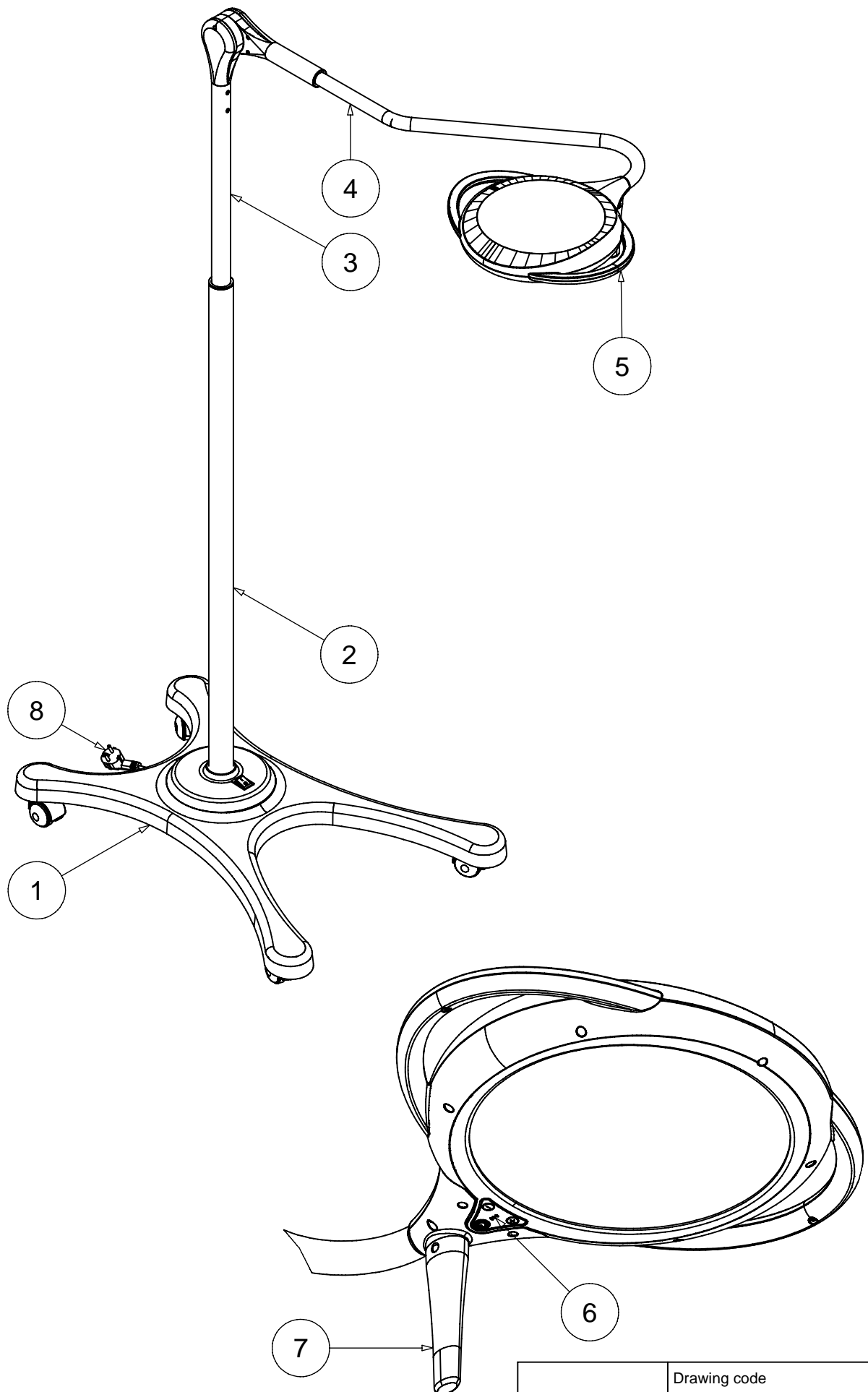
		Drawing code
Rev.	Data	063

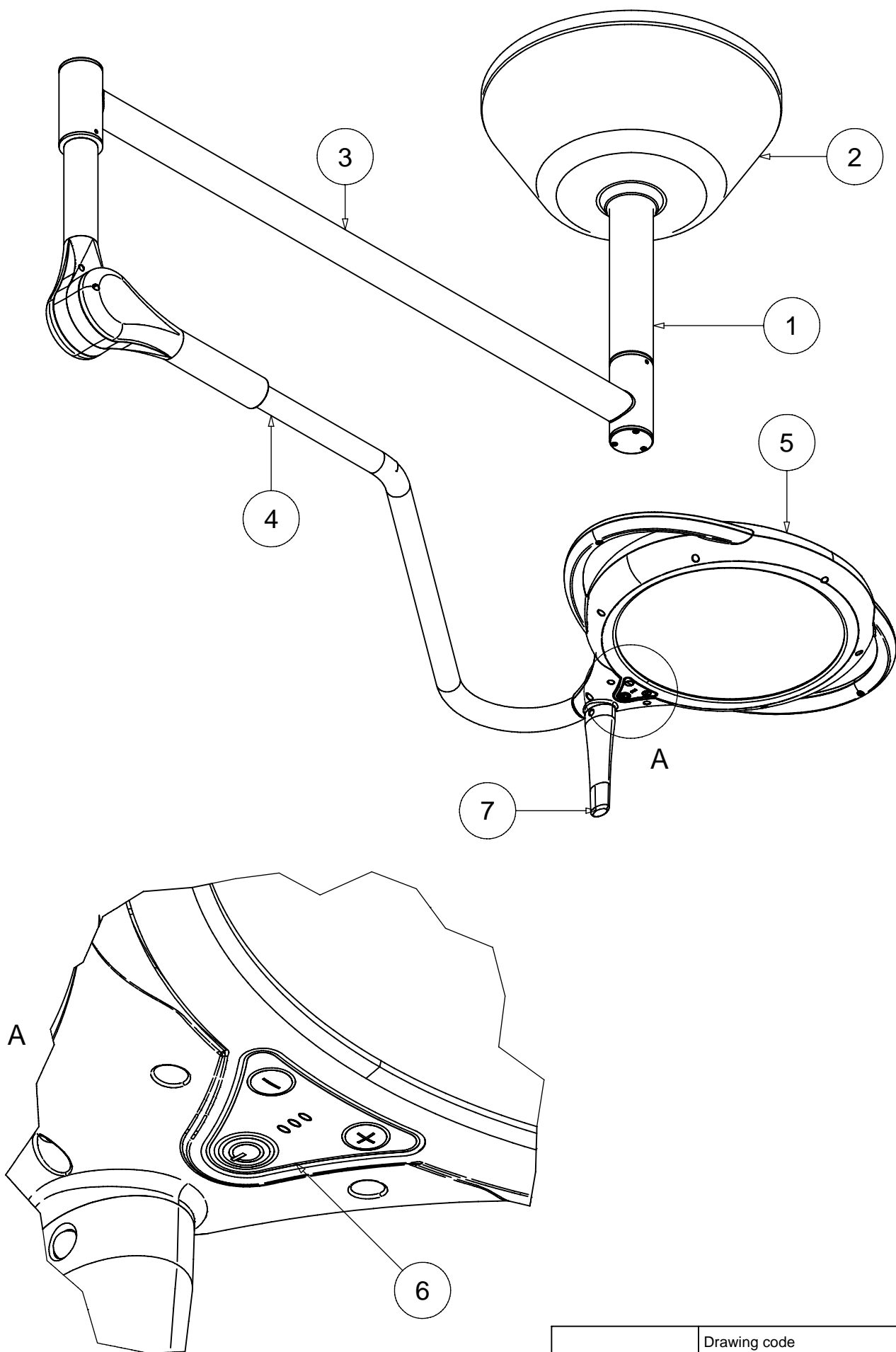		Drawing code
Rev.	Data	064

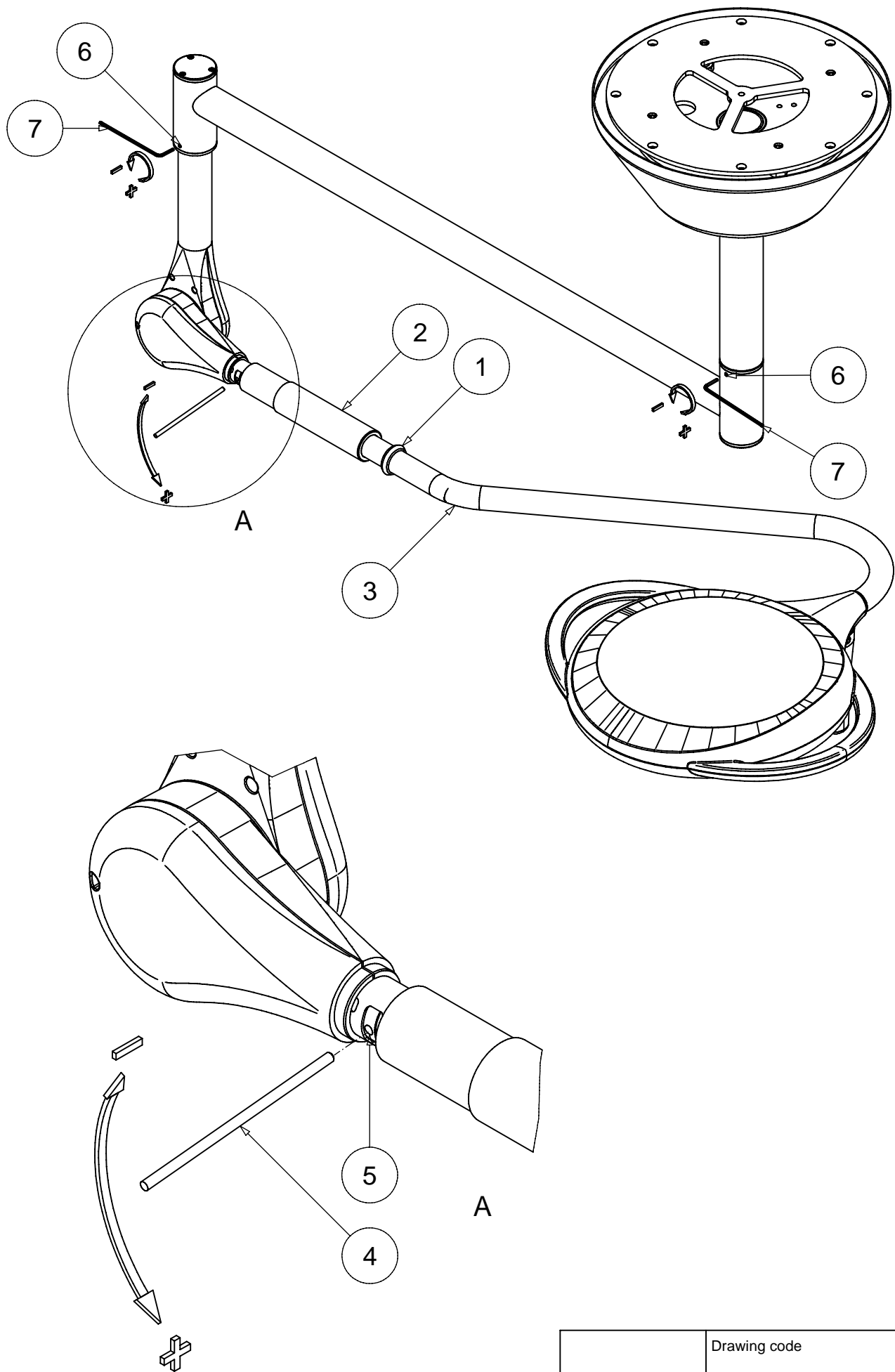		Drawing code	
Rev.	Data	065	

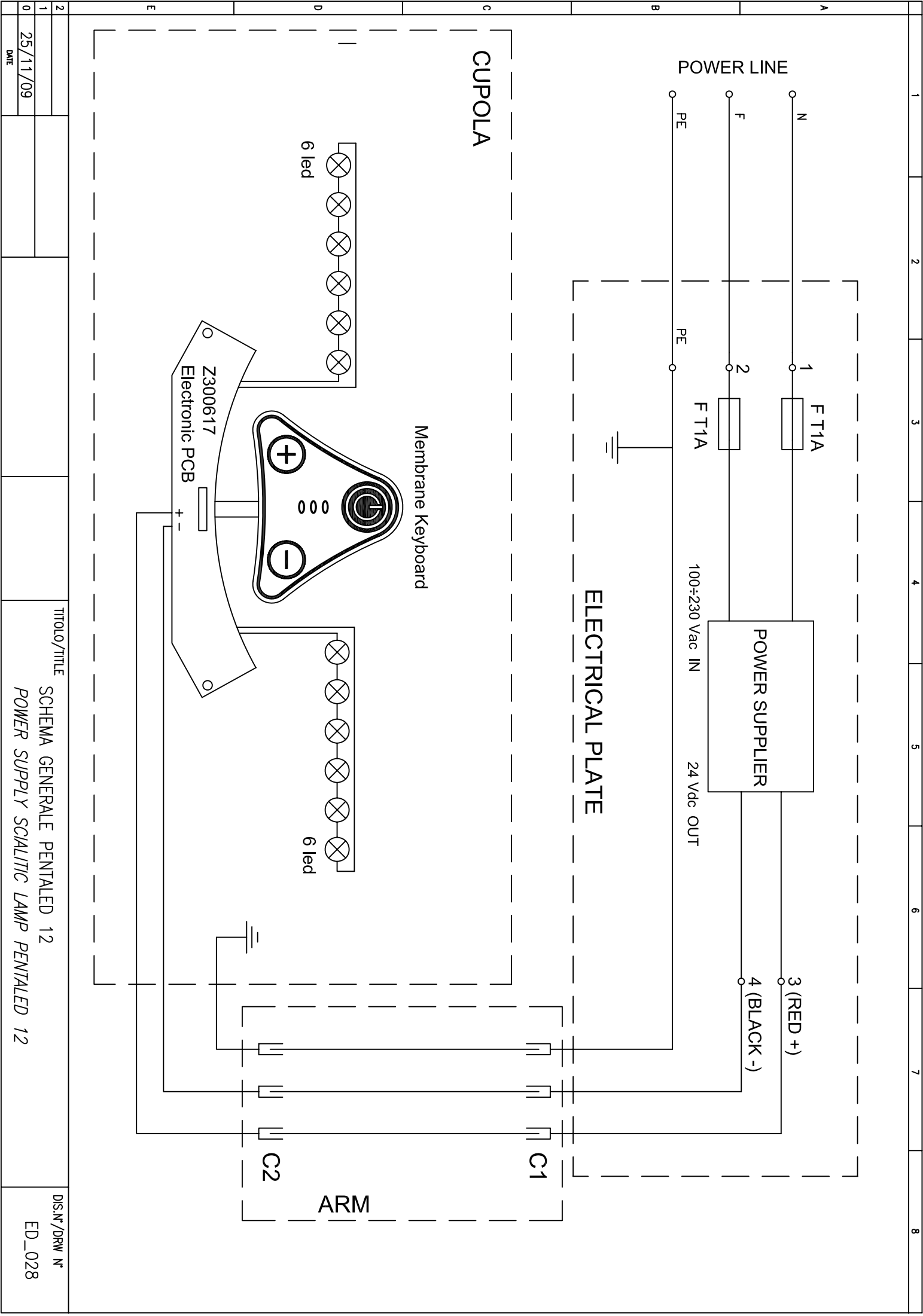2							
1							
0	25/11/09						
	DATE						